




MEDICAL SECTOR

ECG DIGITAL ELECTROCARDIOGRAPH



Kalstein®

A different accompaniment, at your service

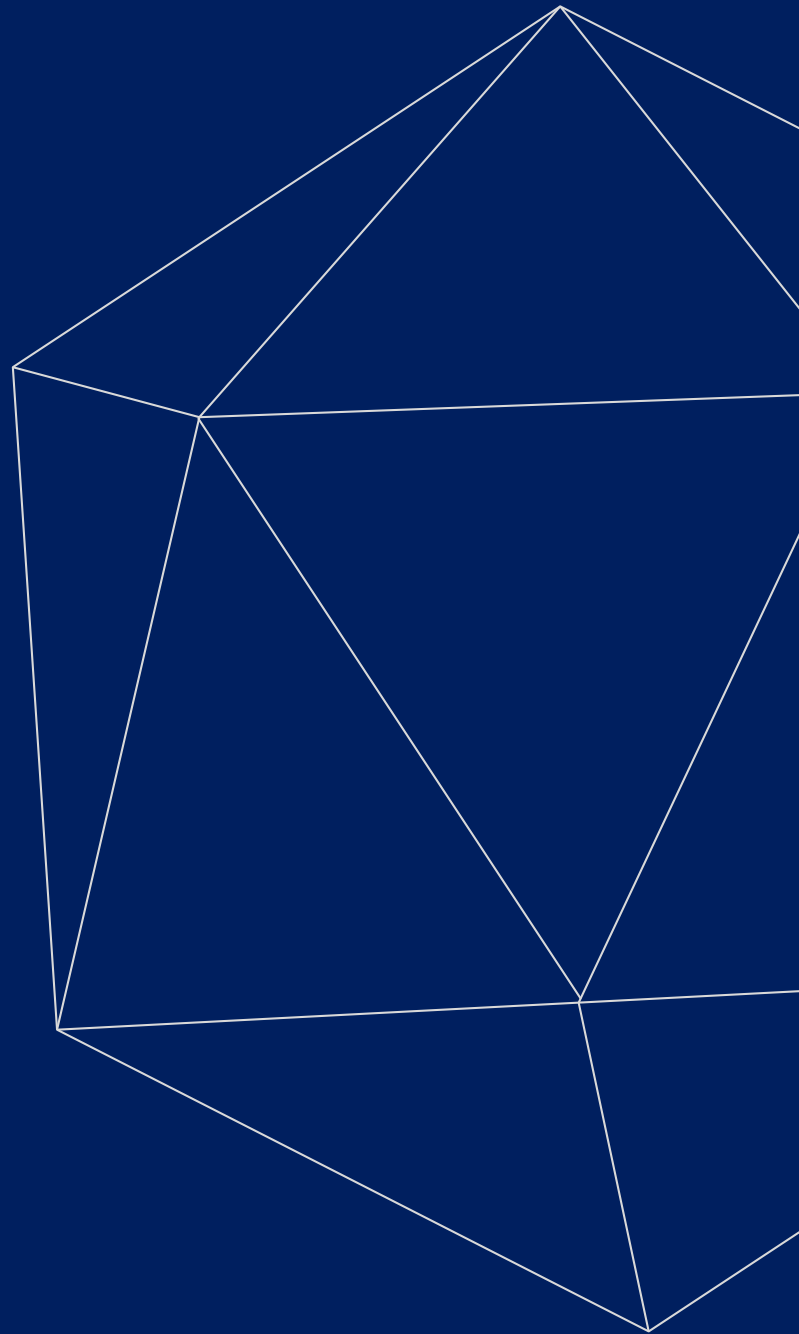



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® ECG Electrocardiograph
12 Channel 12 Lead 6 Channel Medical ECG Machine
Model YR05159



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® ECG Electrocardiograph 12 Channel 12 Lead 6 Channel Medical ECG Machine

Model YR05159



Product Description:

The 12-channel ECG machine is an electrocardiograph that can collect 12-lead ECG signals simultaneously and print ECG waveforms with thermal printing system. The EKG is commonly used to detect abnormal heart rhythms and to investigate the cause of chest pains. An EKG records the electrical activity of the heart. The heart produces small electrical impulses that travel through the heart muscle to cause the heart to contract. These impulses can be detected by the ECG machine.



Product Description:

Model	YR05159
Elements	specification
Characteristics	7 " screen,
	touch screen
	Supports Ecg analysis capabilities
Guide	12 standard conductors
prospect acquisition	synchronously 12 leads
Input circuit	Floating; Defibrillator protection circuit
Input impedance	$\geq 50M\Omega$
Input circuit current	$\leq 0.05\mu A$
Record mode	Automatic: 3CH \times 4 + 1R, 3CH \times 4,
	3CH \times 2 + 2CH \times 3, 3CH \times 2 + 2CH \times 3 + 1R, 6CH \times 2;
	Manual: 3CH, 2CH, 3CH + 1R, 2CH + 1R;
	Rhythm: any selectable lead.
Filter	EMG filter: 25Hz / 30Hz / 40Hz / 75Hz / 100Hz / 150Hz
	DFT filter: 0.05Hz / 0.15Hz AC filter: 50Hz / 60Hz



#Let'sgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® ECG Electrocardiograph 12 Channel 12 Lead 6 Channel Medical ECG Machine

Model YR05159



CMRR	> 100dB
Patient current leak	<10 μ A
Input circuit current	<0.05 μ A
Frequent response	0.05Hz ~ 150Hz (-3dB)
Sensitivity	2.5mm / mV, 5mm / mV, 10mm / mV, 20mm / mV
Anti-baseline drift	Automatic
Constant time	\geq 3.2s
Noise level	<15 μ Vp-p
Paper speed	12.5mm / s, 25mm / s, 50mm / s
Recording mode	Thermal printing system
Paper specification	Sun-7031: 80mm x 20m paper roll Sun-7062: 120mm x 20m paper roll
LCD screen	7 " graphic LCD
Safety rating	IEC60601-1 class I, type CF
Power supply	AC: 100 ~ 240V, 50 / 60Hz, 30VA ~ 100VA
Fuse	AC220V \pm 10% 2 * Φ 5 \times 20mm, T2A / 250V AC time delay



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale

<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

