




LABORATORY SECTOR

ULTRASONIC CLEANER



Kalstein®

A different accompaniment, at your service

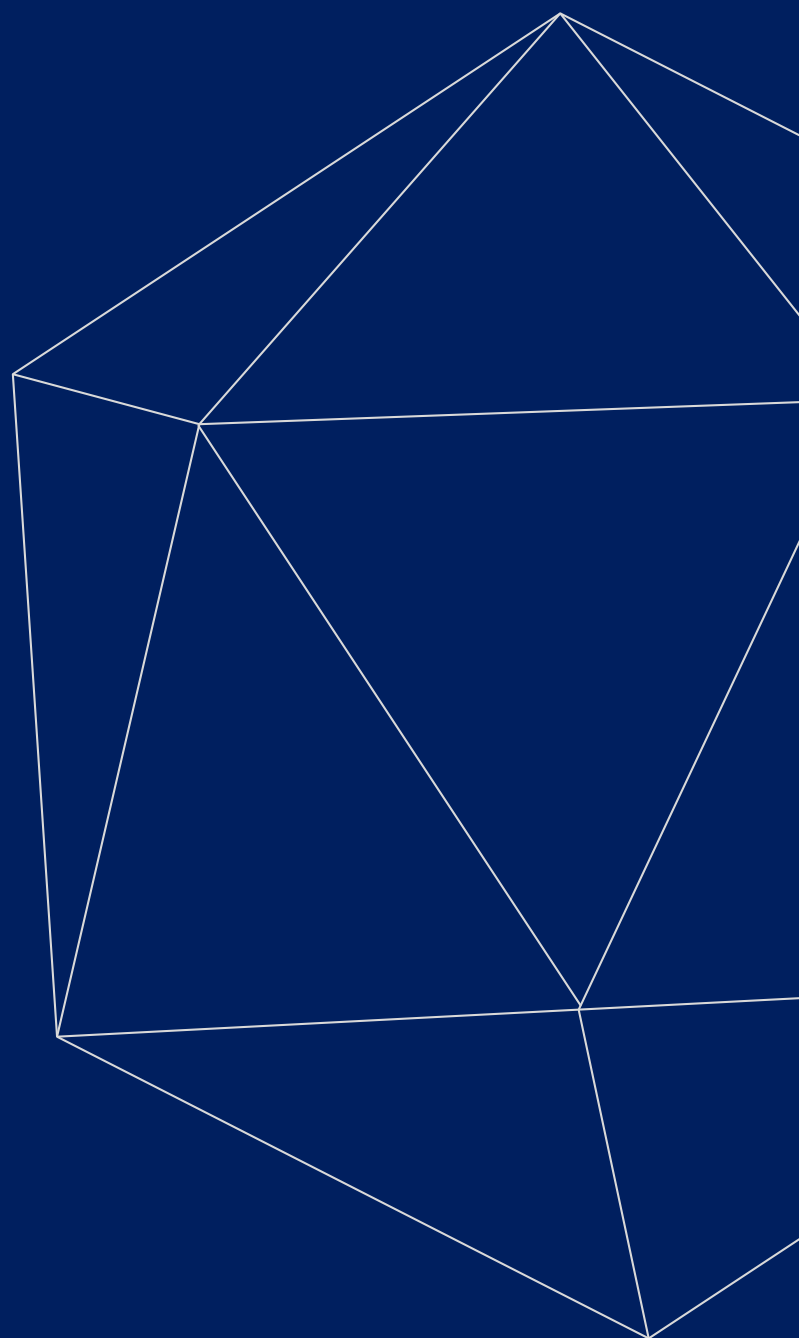



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Ultrasonic Cleaner
2 Liter Ultrasonic Cleaner For Jewelry
Model YR05209



#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S.




KALSTEIN® Ultrasonic Cleaner 2 Liter Ultrasonic Cleaner For Jewelry Model YR05209




Product Description:

YR05209 as a ultimate multi-purpose household cleaner. It's digital timing and heating control, super easy to operate. Clean anywhere you go with this powerful cleaner. As we all know, cleaning jewelry regularly is important and our cleaners can help you out. We also have CE,FCC,ROHS,PSE certificates.

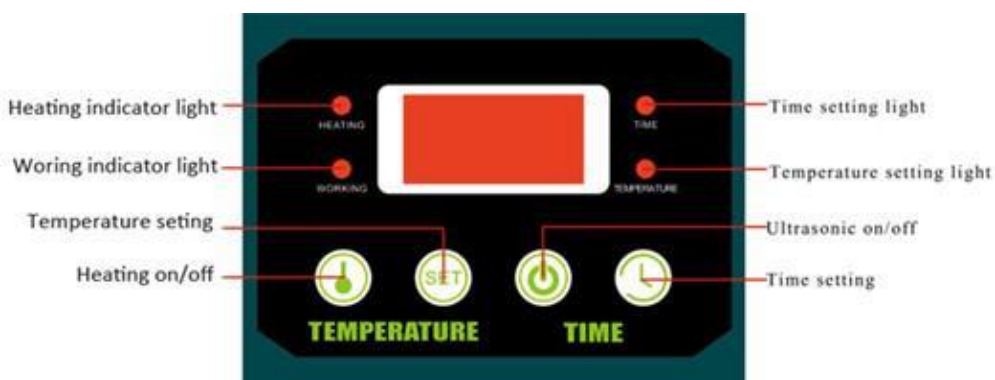

Introduction of Ultrasonic Cleaner:

The high-frequency vibration signal from the ultrasonic cleaner of the digital ultrasonic cleaner in the laboratory is converted into high-frequency mechanical vibration by the transducer, so that it can be spread in the cleaning liquid, and tens of thousands of tiny bubbles will be generated due to the flow of the liquid. The negative pressure area formed in the process of ultrasonic propagation will form and grow, and will close rapidly in the positive pressure area. This process can produce instantaneous high pressure of more than 1000 atmospheric pressure, which is continuous and is called "cavitation effect". It is like a series of explosions, which constantly impact the surface and cracks of objects, so as to achieve a comprehensive cleaning effect..

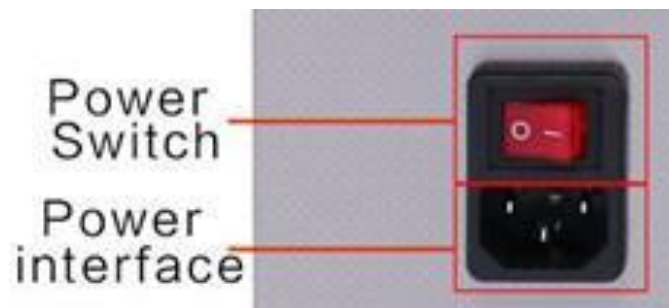

Name	Ultrasonic cleaner	Voltage	<ul style="list-style-type: none"> • AC110-120V 60Hz • AC220-240V 50Hz
Model	YR05209	Power	70W
Packing size	242*227*275mm	Ultrasonic frequency	40KHz
Internal groove size	150*135*100mm	Time setting	10-600S
Tank capacity	2L	Net weight	2.7KG



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

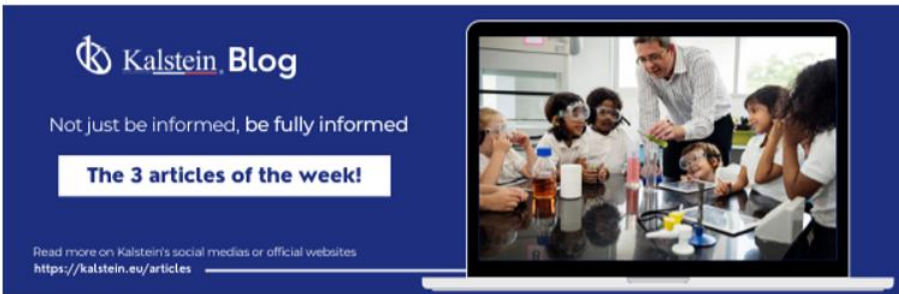

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

