




LABORATORY SECTOR

PH METERS



Kalstein®

A different accompaniment, at your service

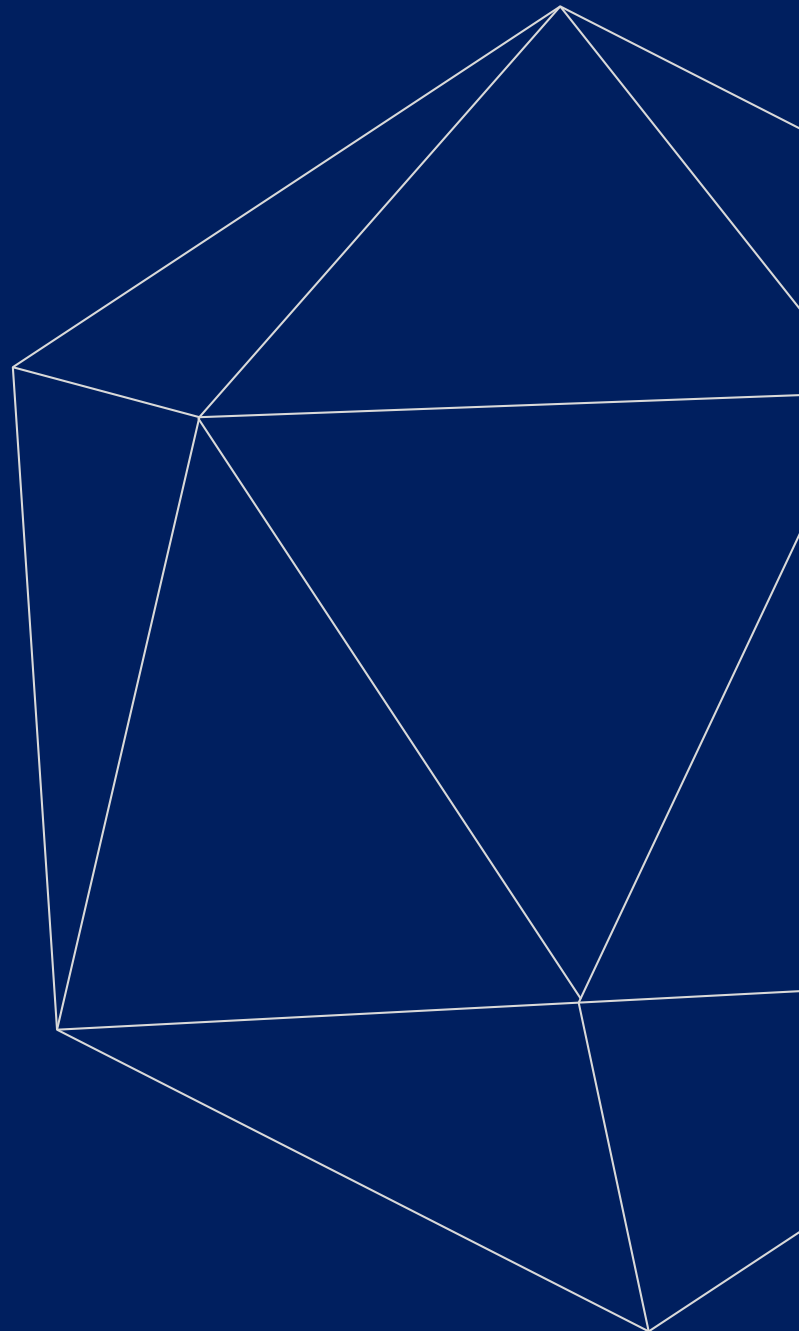



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® PH Meters
Benchtop Dissolved Oxygen Meter
Model YR01834



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® PH Meters
Benchtop Dissolved Oxygen Meter
Model YR01834



Product Description:

Economical benchtop dissolved oxygen meter, 1 or 2 points calibration, setup menu contains 6 options. The meter is suitable for laboratory applications, accuracy: 0.5mg/L.



Features:

- ✓ Economical benchtop dissolved oxygen meter is equipped with a backlit LCD display.
- ✓ 1 or 2 points calibration using the air-saturated water or zero oxygen solution.
- ✓ Salinity and barometric pressure compensations eliminate the measurement error.
- ✓ Automatic temperature compensation ensures accurate readings over the entire range.
- ✓ Auto-Read function senses and locks the measurement endpoint.
- ✓ Setup menu allows setting the number of calibration points, concentration unit, temperature unit, auto-power off, etc.
- ✓ Reset function automatically resumes all settings back to the factory defaults.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® PH Meters

Benchtop Dissolved Oxygen Meter

Model YR01834



Technical Description:

Model	YR01834
Dissolved Oxygen	
Range	0.0 to 20.0mg/L
Resolution	0.1mg/L
Accuracy	±0.5mg/L
% Saturation of Oxygen	
Range	0.0 to 200.0%
Resolution	0.1%
Accuracy	±2.0%
General Specifications	
Calibration Points	1 or 2 points
Temperature Compensation	0 to 40°C/32 to 104°F, automatic
Pressure Correction	60.0 to 112.5kPa/450 to 850mmHg, manual
Salinity Correction	0.0 to 35.0g/L, manual
Sensor Type	DO100 polarographic probe
Connector	6-pin nimi-DIN
Display	Custom LCD (120×60mm)
Power	DC9V, using AC adapter, 220V/50Hz
Dimensions	210(L)×205(W)×75(H)mm
Weight	1.5kg
Order Code	Accessories
YR01834	DO100 dissolved oxygen electrode
	DO electrolyte solution
	DO membrane cap
	DC9V power adapter



#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale

<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89

Cel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

