

LABORATORY LINES CENTRIFUGES







All rights reserved ® KALSTEIN France S. A. S., Optimum Business Center 450 Rue Baden Powell, 34000 Montpellier, France. Tlf: +33 467158849 / +33 680760710/ +33 663810023 https://kalstein.eu KALSTEIN FRANCE, S. A. S





our services Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings In any part of the world, receive your

induction or training from our specialized team of engineers



Quick Response Our work team is always available to response all your consults or questions, in order to support you in any situation.





#Letsgivemore 💚

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or.air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.







Oil Test Centrifuge with LCD Screen Model YR138

Almost every research lab in the world has a centrifuge. They are a necessary equipment to any investigation. That has to do with the fact that the separation process is fundamental when it comes to testing samples and performing important diagnosis. A test centrifuge plays a very important role in many researches in any field.

The centrifuge process happens when the samples start to spin at very high speeds and the components get under the centrifugal force, which causes more dense particles to migrate away depending on the density, viscosity and other qualities.

Oil test centrifuge

This type of centrifuge is designed for testing oils and determines the water and deposit in crude oil. These features make this YR0138 model the ideal equipment to determine moisture in oil exploration and research institutes. This oil test centrifuge is also very common to determine the precipitation number, trace sediments, and insolubility in used lubricating oils.

Oil test centrifuge machine

This device runs with a microcomputer control and it has an AC variable-conversion motor, which delivers high torque, it offers a stable operation, makes very minimum noise and it has a high-speed accuracy. In addition, it has a touching panel with programmable operation, which can be set according to your demands and automatically stored. This particular feature is very important in case of a power failure because the device saves the configuration you want. Another great feature is that this oil test centrifuge comes with an LED screen display with a user-friendly interface with simple and convenient operations. Moreover, this model has great important features such as:

- ✓ Electronic lock.
- ✓ Self-locking device.
- ✓ Over-speed safety device.
- ✓ Over-temperature safety device.
- ✓ Imbalance safety device.
- ✓ Automatic alarm device.

This oil-test centrifuge machine has a real-time conversation of speed and RCF is very convenient for any operation. Finally, this model has 10 acceleration and deceleration files for you to choose.

About this model's speed, you should know that it has a rotor that runs 4000 rpm with a max RCF of 4690g. This test centrifuge's capacity is 100ml / 200mlx4. It has a French Tecumseh Compressor Unit, R404a. In addition, the power supply is AC220V 50Hz 2.5A and the external dimensions are (LxWxH) 650mm x 710 mm x 1100 mm with a net weight of 140 Kg.



CE

Technical Performance

- ✓ Microcomputer control, electronic lock and AC variableconversion motor which delivers high torque, stable operation, low noise, high speed accuracy.
- Touch panel, programmable operation, operation data can be set according to demands and stored automatically.
- ✓ LED display, user-friendly interface, simple and convenient operations.
- ✓ The real-time conversation of Speed and RCF is convenient for operation.
- ✓ It has self-locking device, over-speed safety device, overtemperature safety device, imbalance safety device and automatic alarm device.
- \checkmark It has 10 files of acceleration and deceleration for your choice.



#Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved (B KALSTEIN France S. A. S., Optimum Business Canter 450 Rue Baden Powell, • 34000 Montpellier, France. Tit: +33 467158849 / +33 680760710/ +33 663810023 https://kalstein.eu KALSTEIN FRANCE, S.A. S



www.kalstein.eu

Oil Test Centrifuge with LCD Screen Model YR138



Specifications:

| Model | |
|----------------------|--|
| Max speed | 4000r/min 5000r/min |
| Max RCF | 4690xg 3850xg |
| Max capacity | 100ml/200mlx4 |
| Speed accuracy | ± 50r/min |
| Time setting | 1min~99min |
| Temperature Setting | -20~80°C |
| Temperature Accuracy | ± 2℃ |
| Compressor Unit | French Tecumseh Compressor Unit, R404a |
| Noise | <65dB(A) |
| Power supply | AC110/220V 50Hz/60Hz 2.5kv AC110/220V 50Hz/60Hz 1.5kv |
| Dimension(LxWxH) | 650mmx710mmx1100mm |
| Packing size(LxWxH) | 740mmx830mmx1300mm |
| Net weight | 140kg |

Rotor:



NO.1 Capacity: 100/200ml×4 Speed: 4000r/min Centrifugal force: 4690×g



#Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved (8) KALSTEIN France S. A. S., Optimum Business Center 450 Rue Baden Powell, • 34000 Montpellier, France. Tit: +33 467158849 / +33 680760710/ +33 663810023 https://kalstein.eu KALSTEIN FRANCE, S.A. S





#LETSGIVEMORE

WITH THE **ADCQUISITION OF A KALSTEIN EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:





https://onetreeplanted.org

Humatem Solidarité et coopération biomédicale https://www.humatem.org/



https://www.maniapurefoundation.org/



ANY DOUBT? Contact us!







PARIS - FRANCE **H E A D Q U A R T E R** 2 Rue Jean Lantier 75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89 **Cel:** +33 (0) 6 80 76 07 10 sales@kalstein.eu


https://kalstein.eu/ https://www.kalstein.fr/

Kalstein Training Videos (Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago

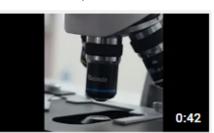

Kalstein's Microscopes 445 views • 1 week ago

