




LABORATORY SECTOR

COLOR LIGHT BOX



Kalstein®

A different accompaniment, at your service

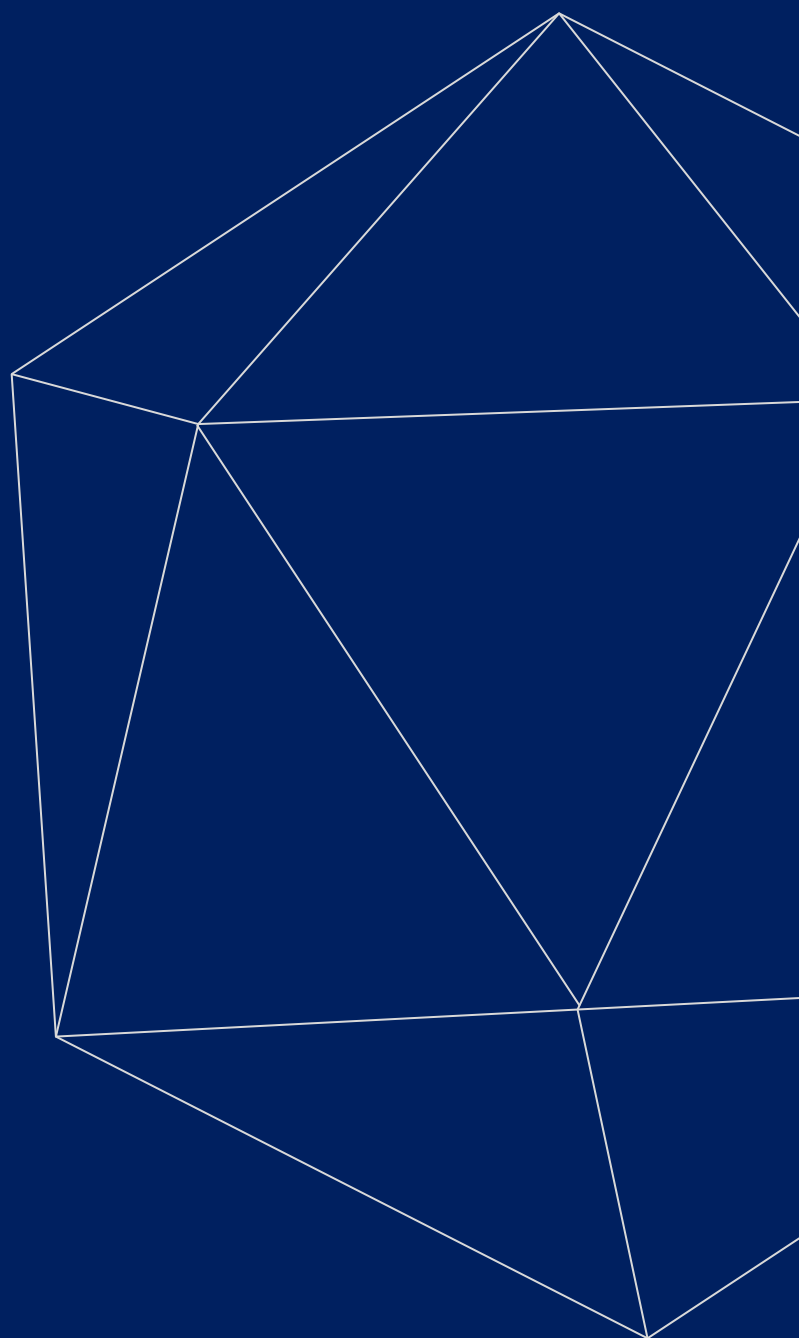



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Color Light Box
Color Proof Station Light Box
Model YR05536 // YR05538



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Color Light Box
Color Proof Station Light Box
Model YR05536 // YR05538



Product Description:

Color Controller Color Viewer color light box is widely applies to printing industry, furnishings and color matching of the printing ink industry etc.. It is a necessity for those industry to control the colors and qualities of products. This Color Viewer color light box supplies D65, D50,U30 light sources. This machine was controlled by microcomputer, which allows the rapid conversion between each light source. The Color Viewer matching box also had the function of display the running time of each light source.



Light Source:

In different condition, the object may present different spectrum with same color. Our vision result of the colors will be differently, thereby it will affect us to judge the color accurately by our naked eyes. While according to our customer's demand, the Color Viewer ,which is the international standard, can provide the right light source you need when you check the colors, such as D65,D50,TL84,UV,U30,CWF etc.. The life-span of lamp depends on the accumulated running time and the turn-on times. When the lamp becomes aging, light instability or the two sides turn black, please replace the lamp timely, otherwise, the checking result will be affected.

Technique Parameter:

Light source: imported fluorescent lamp

- ✓ 36W D65 *4
- ✓ 36W D50 *4
- ✓ 36W U30 *4

Technical Parameters:

Model	YR05536	YR05537	YR05538
Tubes	D65 6500K 36W	D65 6500K 36W	D65 6500K 36W
		D50 5000K 36W	D50 5000K 36W
			U30 3000K 36W
Power	AC220V 50HZ	AC220V 50HZ	AC220V 50HZ
Wortable Size (cm)	135*90	135*90	135*90
Size (cm)	135*99*217	135*99*217	135*99*217



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:

 **Fundación
Jacinto Convit**
<https://www.jacintoconvit.org/>

 **ONE
TREE
PLANTED**

<https://onetreepanted.org>

 **Humatem**
Solidarité et coopération
biomédicale
<https://www.humatem.org/>


**MANIAPURE
FOUNDATION USA**
<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

