




DENTAL SECTOR

DENTAL UNIT



Kalstein®

A different accompaniment, at your service

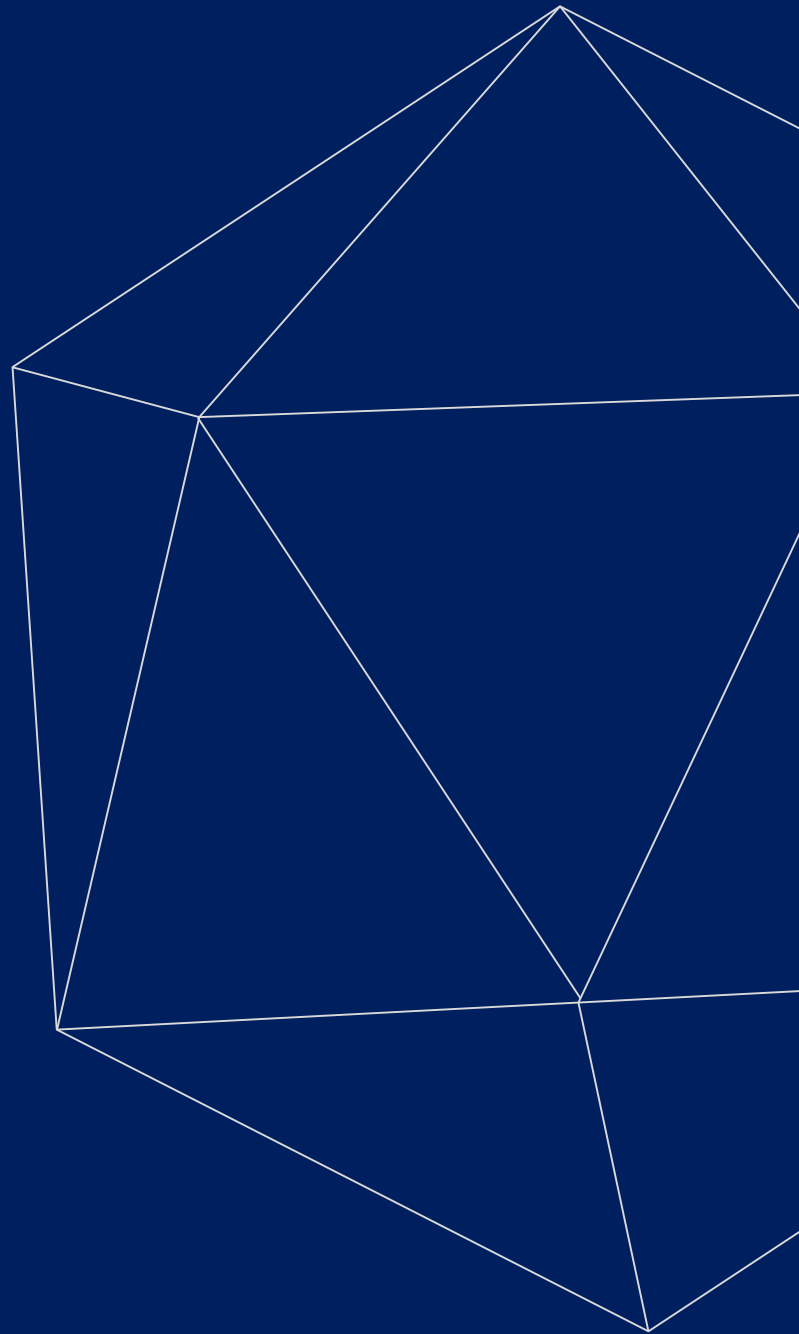



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Dental Units
Disinfection Dental Chair Unit
Model YR06011



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Dental Units
Disinfection Dental Chair Unit
Model YR06011



Product Description:

YR06011 disinfection dental chair, safe and intelligent, one-key disinfection, empowers waterway safety, is All-round coverage, headpiece tube, three-purpose gun, scaled tube, mouthwash pipe and water inlet, completely decompose the biofilm tissue in the pipe. And at the same time eliminate the bacteria and viruses carried by the biofilm, ensuring that the water quality always meets the regulatory requirements.

- ✓ One-button automatic flushing: before the start of the daily diagnosis and treatment, rinse the water outlet of the diagnosis and treatment water for at least 30 seconds, after each diagnosis and treatment, rinse the waterway connected to the oral equipment for at least 30 seconds, which meets the requirements for daily maintenance of the dental unit waterway.
- ✓ Anti-suction device: prevent the water or other media in the instrument from returning to the dental chair price pipeline due to reverse flow.
- ✓ Drainage system: whether it is automatic pipe disinfection or one-button automatic pipe flushing, the position of the chair will automatically rise to the highest position, completely draining the accumulated water in the pipe, preventing the growth of bacteria and odors after the dirt remains.
- ✓ The outer tube of the instrument tube: the pipeline connecting the oral instruments is in direct contact with the ground to prevent the outer tube from being polluted by friction with the floor of the clinic.
- ✓ Highly intelligent disinfection program: one-key disinfection, accurate and efficient.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Dental Units
Disinfection Dental Chair Unit
Model YR06011



- ✓ The medium in the pipe is rinsed with clean water for 2 minutes to start disinfection, which will achieve the highest effect of pipe disinfection.
- ✓ Automatically inject disinfectant water into all oral water pipes, and automatically rinse, saving time and effort.
- ✓ Sterilized water is soaked in the pipeline for 19 minutes and 30 seconds to completely decompose the biofilm tissue in the pipeline.
- ✓ Automatically flush the residual disinfectant in the pipeline for 2 minutes ,and 30 seconds to ensure the safety of the water used for diagnosis and treatment on the human body and the life of the dental chair unit.

Comfortable Dental Treatment Experience:

- ✓ Heavy-duty hydro formed backrest and seat frame, dental chair price gives patients the first sense of security.
- ✓ The ergonomic chair cushion keeps the patient in the most comfortable position at all times. This ensures a more focused and efficient operation by the dentist.
- ✓ Dental chair unit S6 is control system of slowly rising or slowly pausing is adopted to eliminate the patient's frustration caused by the running of the chair.

Standard :

Model	YR06011
3-way syringe(cold/hot)	2sets
Powerful suction device	1set
saliva ejector	1set
Dental light	1set
Warm water system	1set
Pure water supply system	1set
Rotatable ceramic spittoon	1set
DC 24V motor driven chair	1set
Multifunctional assisant device	1set
LED x-ray viewer	1set



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





#LETSGIVEMORE

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**


YOU MAKE A CONTRIBUTION TO:



<https://onetreepanted.org>



Articles Kalstein
(Social Networks and Official Websites)

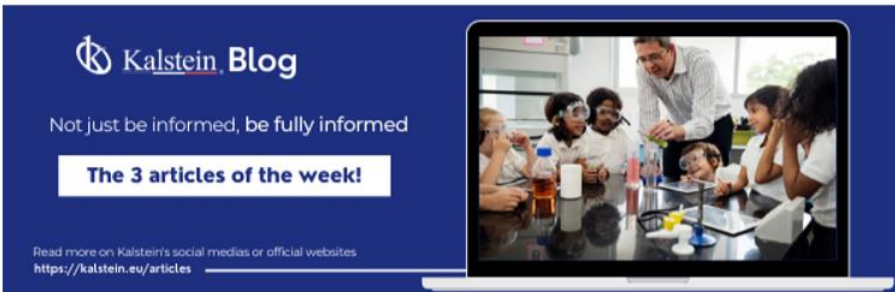

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
microscope is an instrument used for observing objects too small to be seen by the naked eye. This is achieved through the use of lenses or a single reflective surface in which light from the specimen is reflected to form an image.
Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
Refrigeration is a process of removing heat from a space or object to maintain it at a controlled environment (refrigerated space), without fluids and substances.
Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

