




LABORATORY LINES

ICE MAKER



Kalstein®

A different accompaniment, at your service

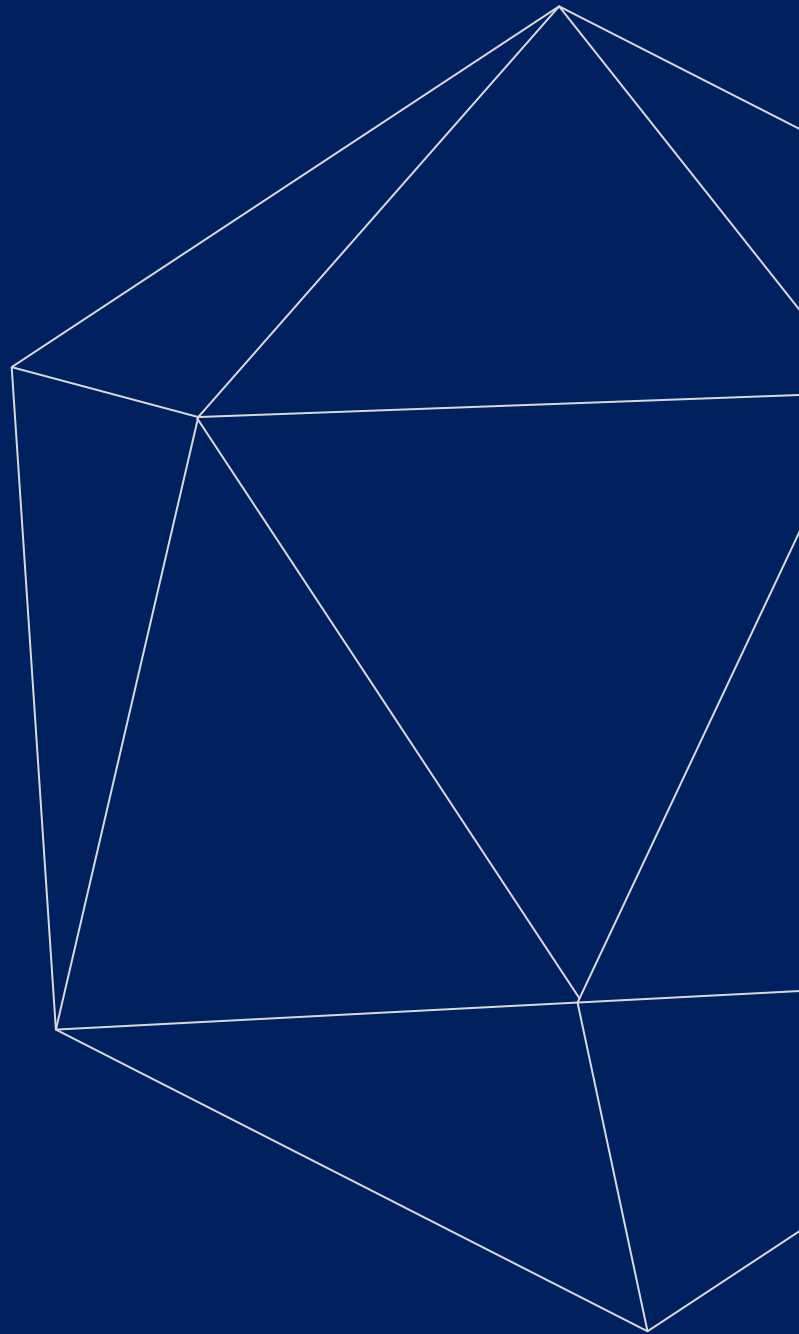



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.





Flake Ice Maker Machine

Series YR05130

Features:

- ✓ High-quality stainless steel is used anti-corrosion and durable. The compact structure is simple and saves space.
- ✓ The cabinet is filled with fluorine-free thermal foam with good insulation effect. The inner lining is food grade plastic. They save energy and are respectful with the environment.
- ✓ World famous brand CFC-free R134a compressors are used. The main parts and components have passed the relevant safety certifications. The product has stable quality and long service life.
- ✓ The ice making process is fully computer controlled, using imported computer chips to make the control reliable and easy to operate.
- ✓ Advanced technologies and imported accessories make low noise, smooth and reliable operation. The vents and fans are deployed to ensure that the geared motor operates reliably.
- ✓ The patented technology of the separate evaporator chip ice chamber line provides high cooling efficiency and large ice making capacity.
- ✓ Unique float type water inlet system ensures no residual water. no need for an ice cleaning process or wasted water, saving water and energy.
- ✓ Long-term performance testing and electrical safety testing will be conducted before shipment, ensuring excellent product performance.
- ✓ Fully automatic operation with full ice screen, water shortage screen, operation screen and malfunction warning screen, etc. For power-off protection. They will stop automatically in case of ice or water shortage and start automatically when the ice is not full or the water supply is restored.

Ice Features:

- ✓ Screw extrusion hob style ice making type makes compact structure and automatic water ice separation. The optimal design of the skate blade makes the ice shape small and practical.
- ✓ Flake ice cools faster than other forms of ice, which can infiltrate a narrow space, with fast cooling speed, good ice bath effect.



#Letsgivingmore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Flake Ice Maker Machine Series YR05130

 Energy-saving and environmental friendly



Parameters per Model

Model	YR05130	YR05131
Ice Making Capacity	20 Kg / 24h	30 Kg / 24h
Ice Storage	10Kg	10Kg
Voltage	220V/50Hz or 110V/60HZ (should pay extra charge)	
Power	280 W	300 W
Cooling Type	Air Cooling	
Refrigerant	R134a	
Water Supply	Tap Water	
Net Weight (Kg)	26	
Gross Weight (Kg)	29	
Product Dimension	330*470*610 mm	
Packagind Dimension	400*550*660 mm	
Stuffing Qty	20'GP	168
	40'GP	348
	40'HQ	464



#Letsgivemore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Flake Ice Maker Machine

Series YR05130

01

Thickened Insulation Layer


Strong insulation bin, ice can be stored for 5-8 hours after ice full shutdown.



02

SUS304

High-quality 304 stainless steel, beautiful and durable.

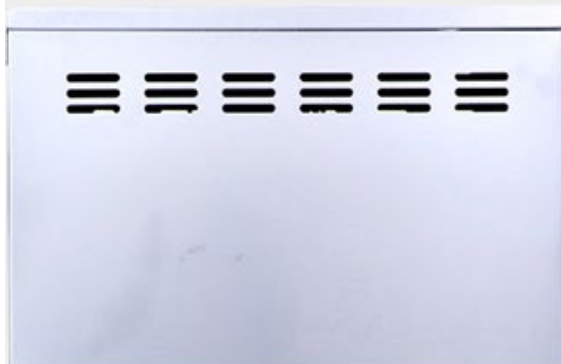

03

Water Filter

Set the filter at the inlet, to prevent large particles into the machine.



04

All-round Cooling

Cooling holes are installed around the machine, to avoid high temperature operation of the compressor.



#Letsgivemore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



 Energy-saving and environmental friendly



Flake Ice Maker Machine

Serie YR05130



#Letsgivingore 


Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Titl: +33 467158849 / +33 680760710 / +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Articles Kalstein
(Social Networks and Official Websites)

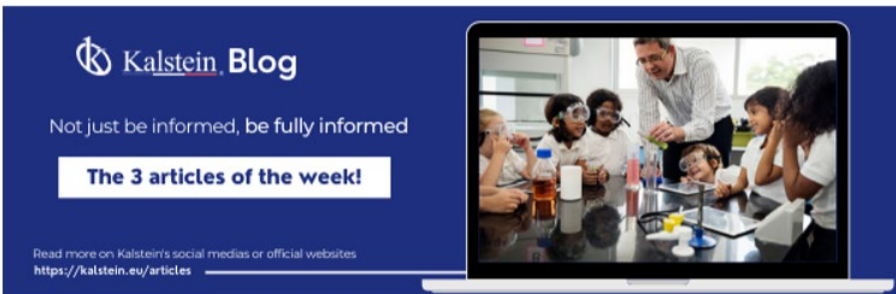

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator in a laboratory is one of the most important pieces of equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances...

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tif: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

