




# LABORATORY SECTOR

MICROSCOPE



## Kalstein®

A different accompaniment, at your service

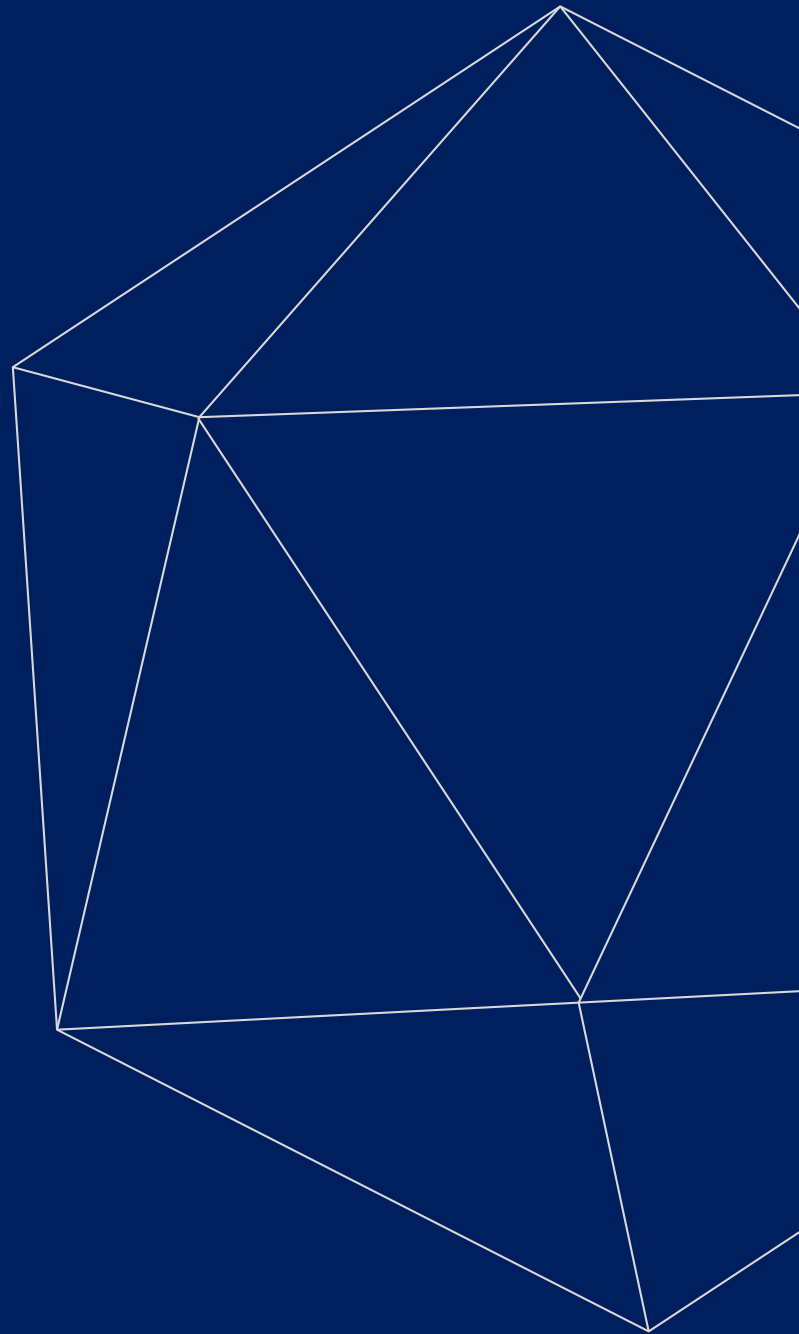



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



**KALSTEIN® Microscope**  
**Laboratory Biological Microscope, Trinocular Achromatic**  
**Model YR05788-AT**



**#Letsgivemore** ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
**KALSTEIN FRANCE, S. A. S**



**KALSTEIN® Microscope**  
**Laboratory Biological Microscope, Trinocular Achromatic**  
**Model YR05788-AT**



**Product Description:**

- ✓ Entry Level Laboratory Microscope With Full Range Upgrade Extension.
- ✓ Can Use 40x LWD Infinity Plan Objective W.D.=3.0m To Avoid Contaminated By Cedar Oil.
- ✓ Creative 100x Water Objective Instead Of 100x Oil Objective For Easy Use Without Oil.
- ✓ High Rigidity Body Structure Design Supply High Steady & Easy Handle.
- ✓ Different Color Available To Customize Your Model!



**Seidentopf Head:**

Inclined 30°, rotatable 360°, interpupillary distance 48-75mm, high eye point eyepiece, which is convenient for people with different sitting height; the integrated design of the eyepiece lens barrel, avoids the loss and contamination of the eyepiece and enhances the integrated performance of the microscope.

**Quadruple Backward:**

4 sets objective available: achromatic, plan, infinity e-plan, infinity plan.

**Double Slide Holder Working Stage:**

Mechanical stage size: 145 mm×140mm, travel: 76mm×52mm, scale: 0.1mm.

**Condenser and Light Source:**

Abbe condenser n.a.1.25, iris diaphragm, rack & pinion up/down, with filter holder 3w-led illumination systems.



**#Letsgivemore** ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
 2 Rue Jean Lantier • 75001 Paris •  
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
 KALSTEIN FRANCE, S. A. S



**KALSTEIN® Microscope**  
**Laboratory Biological Microscope, Trinocular Achromatic**  
**Model YR05788-AT**



**Technical Description:**

<b>Model</b>	YR05788-AT
<b>Power Supply Type</b>	DC USB Adapter
<b>Transmit Light</b>	3W LED
<b>Light Source Type</b>	LED
<b>Condenser N.A.</b>	1.25
<b>Condenser Type</b>	Abbe
<b>Condenser Features</b>	Rack & Pinion↓↑, Filter Holder, Iris Diaphragm
<b>Focusing Knob</b>	Low Position Focusing Knob For Easy Operation
<b>Focusing Scale</b>	0.002mm
<b>Focusing Range</b>	26mm
<b>Focusing Type</b>	Coaxial Coarse & Fine
<b>Focusing Features</b>	Safety Stop Screw, Tension Adjustable
<b>Working Stage Moving</b>	76x52mm
<b>Working Stage Size</b>	145x140mm
<b>Working Stage Type</b>	2 Layers Mechanical
<b>Working Stage Features</b>	Double Slide Holder, Size ≥140x140mm
<b>Objective</b>	Achromatic Objective 4x,10x,40x,100x
<b>Objective Type</b>	Achromatic
<b>Nosepiece Holes</b>	4 Holes
<b>Nosepiece Features</b>	Inward
<b>Eyepiece</b>	WF10x/18mm, 2 pcs
<b>Eyepiece Diameter</b>	23.2mm
<b>Eyepiece Features</b>	Diopter Adjustable
<b>Head Light Path</b>	E80:P20
<b>Head Diopter</b>	Left Adjustable
<b>Head Interpupillary</b>	48~75mm
<b>Head Rotation</b>	Head Rotate 360°
<b>Head Inclination</b>	30°
<b>Head Type</b>	Seidentopf Trinocular
<b>Standard Magnification</b>	40x-1000x
<b>Optical System</b>	Achromatic
<b>Upgrade Function</b>	BF+PH+DF+PL+FL
<b>Standard Function</b>	BF



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
 2 Rue Jean Lantier • 75001 Paris •  
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
**KALSTEIN FRANCE, S. A. S**







 **Kalstein**<sup>®</sup>  
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE  
ACQUISITION  
OF A KALSTEIN  
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación  
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



**Humatem**

Solidarité et coopération  
biomédicale

<https://www.humatem.org/>




**MANIAPURE  
FOUNDATION USA**

<https://www.maniapurefoundation.org/>



**Articles Kalstein**  
(Social Networks and Official Websites)

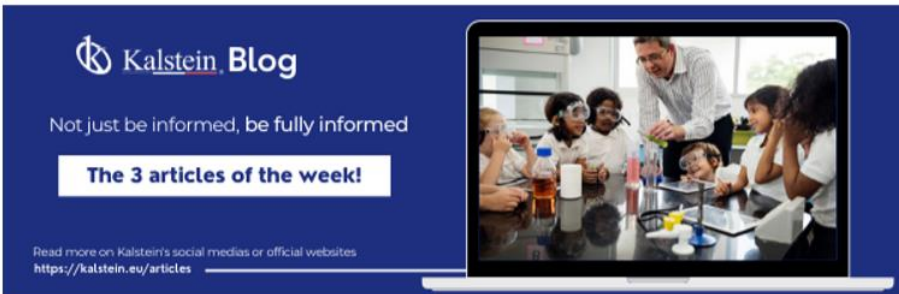

**Kalstein Blog**  
Not just be informed, be fully informed

**The 3 articles of the week!**

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW

**How does the climate affect your microscope?**

**Kalstein Blog**

**Article of the week:**  
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW

**What is the importance of the laboratory refrigerator?**

**Kalstein Blog**

**Article of the week:**  
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE  
**HEADQUARTER**

2 Rue Jean Lantier  
75001 Paris - France

**Tlf:** +33 (0) 01 78 95 87 89

**Cel:** +33 (0) 6 80 76 07 10

[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

