




# LABORATORY LINES

## ICE MACHINE



**Kalstein**®

A different accompaniment, at your service

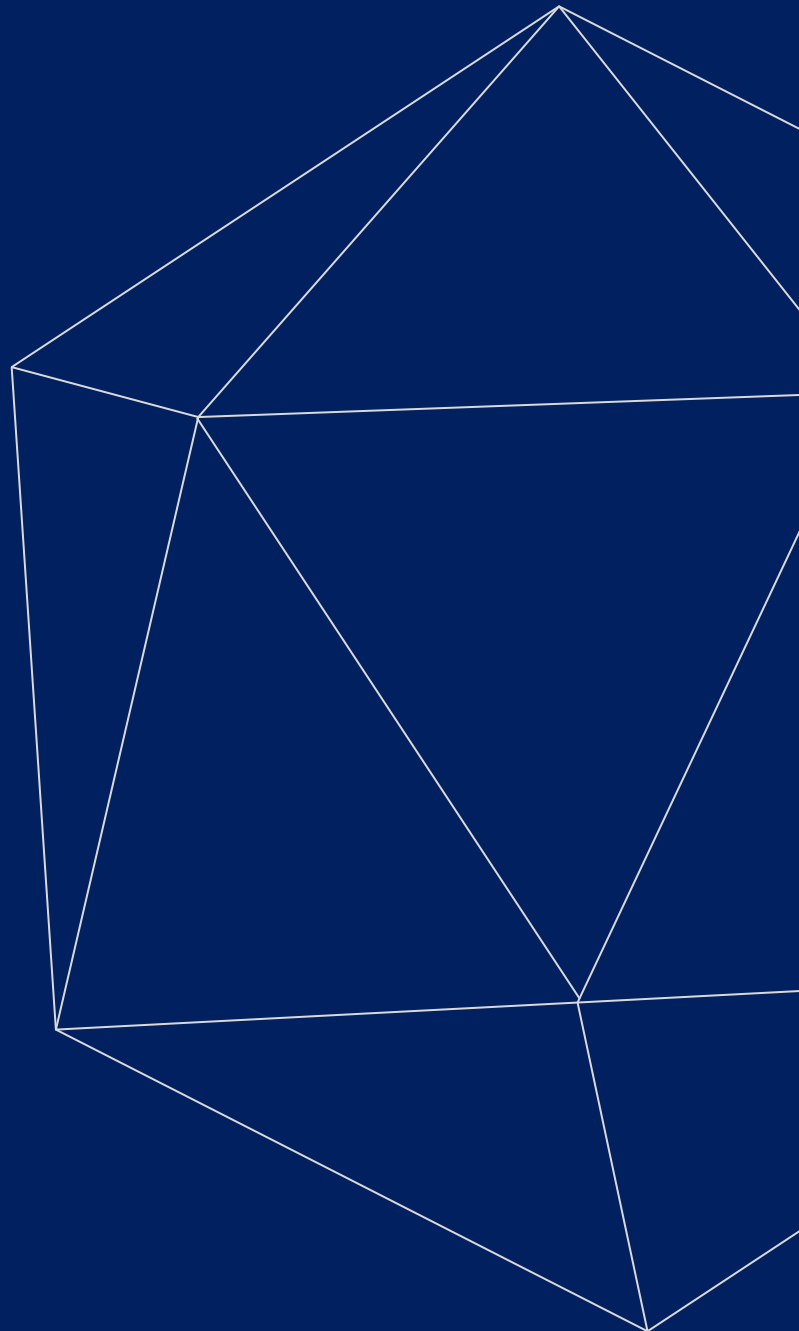



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



## KALSTEIN® Snowflake Ice Machine

Series YR05137 / YR05138 / YR05139



### Features:

- ✓ Use high-quality stainless steel shell, corrosion resistant durable, All-in-one structure, compact and simple, save space.
- ✓ Machine insulation layer is non-fluorine foaming good insulation effect, inner bile is non-fluorin antiseptic type, energy saving and environment protection.
- ✓ Use high-quality and efficient fluorine-free compressor, the main components have relevant safety certification, product safety and security.
- ✓ Ice-making process using all-computer program control, import computer chips, reliable control smooth operation.
- ✓ Use foreign imports of gearboxes and motors, low noise, smooth and reliable operation. The top of the ice machine is equipped with cooling holes and fan to ensure that the reducer motor can also operate reliably in high temperatures and harsh conditions.
- ✓ Use patented technology of compartment ice evaporator, high cooling efficiency, large amount of ice production.



### Ice Features

- ✓ Spiral roller extrusion ice type, compact structure, realistic ice, water automatic separation. The optimized design of the ice blade is to make the ice shape small and practical.
- ✓ Unique water tank floating water intake system, to ensure no residual water, no de-icing process, no water loss, water saving and energy saving.
- ✓ There is ice full display, water shortage display, too cold protection display, fault warning display and other protective shutdown functions. Ice machine will stop working when ice full or water shortage, when water inlet will automatically power on, with automatic memory recovery function.
- ✓ The ice shape is irregular snowflake crushing ice, ice-shaped small, can penetrate the narrow gap, cooling speed, ice bath effect is good, designed for the laboratory.
- ✓ The front is equipped with power switch and function indicators, thoughtful detailed operating instructions, the use of intuitive and convenient, all safety indicators have been tested by electrical performance, safety and rest assured.



#Letsgivemore ❤️

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Tél: +33 467158849 / +33 680760710 / +33 663810023  
<https://kalstein.eu>  
KALSTEIN FRANCE, S.A. S



## KALSTEIN® Snowflake Ice Machine

Series YR05137 / YR05138 / YR05139



### Parameters per Model

Model	YR05137	YR05138	YR05139
Ice Making Capacity	100 kg/24h	120 kg/24h	150 kg/24h
Ice Storage Capacity	40 kg		
Ice Shape	Irregular Small Granular		
Power/ Voltage	520w/220v	550w/220v	680w/220v
Product Dimension	500×611×940mm		
Packing Details	570x680x1000		
N.W / G.W	58 KG / 63 KG		



**#Letsgivingmore** ❤️

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Tél: +33 467158849 / +33 680760710/ +33 663810023  
<https://kalstein.eu>  
KALSTEIN FRANCE, S.A. S







 **Kalstein**<sup>®</sup>  
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE  
ACQUISITION  
OF A KALSTEIN  
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:

 **Fundación  
Jacinto Convit**  
<https://www.jacintoconvit.org/>


  
<https://onetreepanted.org>

 **Humatem**  
Solidarité et coopération  
biomédicale  
<https://www.humatem.org/>

  
**MANIAPURE  
FOUNDATION USA**  
<https://www.maniapurefoundation.org/>



**Articles Kalstein**  
(Social Networks and Official Websites)

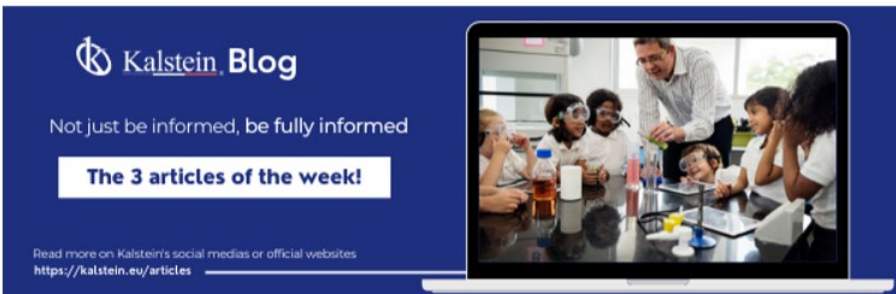

**Kalstein Blog**  
Not just be informed, be fully informed  
**The 3 articles of the week!**  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW  
**How does the climate affect your microscope?**  
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

**Kalstein Blog**  
**Article of the week:**  
How does the climate affect your microscope?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW  
**What is the importance of the laboratory refrigerator?**  
A refrigerator is a device that is one of the most common pieces of laboratory equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances...

**Kalstein Blog**  
**Article of the week:**  
What is the importance of the laboratory refrigerator?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



**Kalstein France - Jacinto Convit Foundation / Alliances**  
23 views • 3 years ago



**Kalstein's Microscopes**  
445 views • 1 week ago



MONTPELLIER - FRANCE  
**HEADQUARTER**  
Optimum Business Center,  
450 Rue Baden Powell,  
34000 Montpellier - France.

**Fax:** +33 (0) 4 67 15 88 49  
**Tif:** +33 (0) 6 80 76 07 10  
[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

**Official Websites for Latin America.**

